

MONOSTENT™


MONOCANALICULAR STENT



THE ONE CHOICE FOR REPAIR OF LACERATED CANALICULI

MonoStent monocanalicular stent is a 50mm length of medical grade silicone tubing with a modified EagleVision FlexPlug™ for immobilization at one end. It comes pre-loaded on a malleable stainless steel stylet making insertion quick and easy. This gives MonoStent enough rigidity to be advanced as far as the nasolacrimal duct if desired.

The opaque FlexPlug allows the MonoStent to be visually monitored. EagleVision designed the unique contouring traction ribs on the tapered shaft of the FlexPlug and the low profile rim to provide superior patient comfort and to minimize extrusion, so that stitching may not be required. In most adult cases MonoStent can be removed in the office.

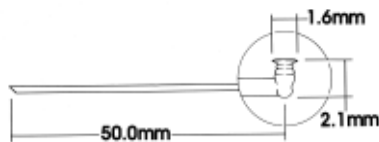

- 50mm length of medical grade silicone tubing
- FlexPlug™ for immobilization at one end
- Preloaded on malleable stylet
- Superior patient comfort

Supplied	Description	REF
1 Stent preloaded on a disposable stylet	MonoStent Monocanalicular Stent	1085

Patents Pending



8500 Wolf Lake Dr., Suite 110 • P. O. Box 34877 • Memphis, TN 38133 USA
 1-800-222-7584 • (901) 380-7000 • Fax (901) 380-7001
www.eaglevis.com • info@eaglevis.com